**WWD700**

**DIGITAL BABY SCALE**  
**PESE-BEBE NUMERIQUE**  
**Dijital BEBEK TARTISI**

**USER MANUAL** EN

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**KULLANIM KILAVUZU** TR

# Contents

|   |          |
|---|----------|
| <b>Technical Properties .....</b>         | <b>3</b> |
| <b>Button Functions .....</b>             | <b>4</b> |
| ON/OFF/TARE Button: .....                 | 4        |
| UNIT/MUSIC Button: .....                  | 4        |
| SET Button:.....                          | 4        |
| <b>Special Operation Functions: .....</b> | <b>4</b> |
| Tare Off:.....                            | 4        |
| Weighing Function:.....                   | 4        |
| Playing Music Function: .....             | 5        |
| Height Measuring Function:.....           | 5        |
| Automatic Screen Backlight: .....         | 5        |
| <b>Fast and Easy Operation .....</b>      | <b>5</b> |
| <b>Cautions .....</b>                     | <b>5</b> |
| <b>Important!.....</b>                    | <b>6</b> |

# Technical Properties

1. *Maximum Capacity : 20kg*
2. *Accuracy Resolution : 5g*
3. *Minimum Weight : 100g*
4. *Units of Measure : kg/lb*
5. *Maximum Tare Weight : <10kg*
6. *Low Battery Indicator: When battery power level is lower than 4,5+/-0,1V, the "Lo" warning sign becomes visible on the screen and the scale automatically powers off.*
7. *Overweight indicator: If Baby's weight is more than 20, 10 kg, the "O-Ld" warning sign becomes visible on the screen.*
8. *Out of Balance Indicator: If the scale does not settle on a smooth, flat surface, the "UnST" sign becomes visible on the screen instead of 0,000kg when the scale powers on.*
9. *Operation Current: On "OFF" State  $\leq 10\mu\text{A}$ , Weighing State  $\leq 6\text{mA}$ , Backlight Active State  $\leq 15\text{mA}$ , Music on State  $\leq 80\text{mA}$*
10. *Supply Voltage: 6V, with 4 \* AA batteries*
11. *Automatic Power Off: If the scale do not use for 60 second, it automatically powers off.*
12. *Accuracy : 10kg+/-10g+1%, 20kg+/-15g+1%, 30kg+/-20g+1%*
13. *Function Buttons: The scale have 3 function buttons, they are ON/OFF/TARE, UNIT/MUSIC and SET*
14. *Backlight on the Screen: It is automatically blue, and not able to change.*
15. *Playing Music Function: There are 20 different songs.*
16. *Ambient Thermometer: 5-50 °C*
17. *Height Measuring Function: 150 cm measuring tape is provided with the device.*
18. *Appropriate Operating Temperature: 10-40°C*
19. *Storage Temperature: -10,50°C*
20. *Humidity :  $\leq 90\%$  Relative Humidity*

# Button Functions

## **ON/OFF/TARE Button:**

- *OFF state: Press this button to power on.*
- *ON state: Press this button to Tare off.*
- *ON State: Hold this button for a few seconds to power off.*

## **UNIT/MUSIC Button:**

- *Press this button to play music while weighing, and press the same button to change music.*
- *Hold this button for a few seconds to change unit.*

## **SET Button:**

- *While weighing, press SET button to see previous weight measurements. It switches back to weighing state in three seconds. The scale stores last 10 measurements.*

# Special Operation Functions:

## **Tare Off:**

1. *While the scale is off, put something on the scale less than 10 kg like blanket, then press ON button, get back the blanket after see 0,000 kg on screen. You see a minus value on the screen after that, you can delete by pressing TARE Button.*
2. *Turn on the scale; put something on the scale less than 10 kg like a blanket when "0,000 kg" is displayed on the screen. Screen shows the weight of blanket. When the value fixed, press TARE Button, this provide to see zero value on the screen. If you get back the blanket, a minus value appears on the screen. If you press TARE Button, you see zero value on the screen.*

## **Weighing Function:**

Put the scale on a horizontal plane. After turned on and saw 0,000 kg on the screen, place your baby on the center of baby scale. After the value is fixed, music begins to play. You can use weighing function again after music stops.

## **Playing Music Function:**

1. *Playing Music; while the scale is on, press the MUSIC Button to play music. You can change music by pressing the same button. There are 20 music options.*
2. *Music ON/OFF; while the scale is OFF, pres and hold SET button and pres ON at the same time. You will see "ON" on the screen, music function is enabled. When you pres the SET key again, you will see "OFF" on the screen and music function is disabled.*

## **Height Measuring Function:**

There is a 150cm measuring tape on the left side of the scale to measure your baby's height. You can rewind measuring tape by pressing the button located under the measuring tape exit hole.

## **Automatic Screen Backlight:**

If there is no operation on the scale, backlight will turn off automatically after 15 seconds. It will turn on automatically when the scale is re-operated.

## **Temperature Switch Between C and F:**

Press button "ON" to start the device (there will be 8888 on screen. Immediately press button "UNIT". You can change the unit of the temeperature sign from C to F ,or from F to C.

## **Fast and Easy Operation**

1. *Open the package and insert the batteries.*
2. *Put the scale on a horizontal surface, and lay a blanket on it. Do not let the blanket touch the ground or extend over the edge of the scale, it might cause inaccurate weighing.*
3. *Press ON/OFF/TARE Button and wait to see 0,000kg on the screen.*
4. *Place your baby on the center of the scale.*
5. *You can see the weight of your baby when the value is fixed and music starts to play (please see playing music function above).*
6. *In case there is no operation the scale turns off in 60 seconds automatically.*

## **Cautions**

1. *Place your baby on the center of the scale for an accurate weighing.*
2. *If there is anything unusual, please remove the battery and wait for 10 minutes. Then replace batteries and try again.*

3. *If “Lo” sign appears on the screen, please change your batteries.  
(Alkaline batteries perform better)*
4. *Keep your scale clean and do not use corrosive liquids and detergents.*
5. *Do not use a load 1, 5 times heavier than maximum capacity of the scale. Use the device on a smooth, horizontal plane and protect it from external shocks.*
6. *If you do not use device frequently, remove the batteries.*
7. *If the scale is used in a long period and frequently, ambient humidity should be  $\leq 90\%$ . (No condensation)*
8. *If batteries have low charge or no charge, please replace batteries to avoid leakage or corrosion on metal contacts.*
9. *If a problem, which cannot be solved by reading the user manual, persists, please remove the batteries and wait for 10 minutes. Then replace the batteries and try again.*

## **Important!**

- *Please make sure total load does not exceed 20 kg. Otherwise you may damage the scale.*
- *Please do not drag the scale. If you want to change its location, please completely lift the device off the ground and place it gently on a new surface.*
- *The scale should be placed on a smooth and flat surface, do not operate the device while carrying it.*
- *Please do not use it with wet hands.*

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### **Üretici:**

BEYON INDUSTRY CO., LTD  
6th Industrial Park, Wusha Village,  
Chang'an Town, DongGuan City, China  
Tel: +86-769-33238965  
Fax:+86-769-89032220-1413

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## Beyon Industry CO.,LTD

ADDRESS: No:3, Zhanhong Street, Dongcheng Road Dongshen town, Zhongshan City, China

TEL NO..+86-769-33238965

FAX NO.+86-769-89032220

### EC Declaration of Conformity

We, Beyon Industry CO.,LTD. No:3, Zhanhong Street, Dongcheng Road Dongshen town, Zhongshan City, China

Declare under our sole responsibility that the product(s):

**Description** : Baby Scale

**Brand Name** : WEEWELL

**Model Number** : WWD700

To which this declaration relates is in conformity with the essential requirements and other relevant requirements of EMC Directive (2014/30/EU). The product is in conformity with the following standards and/or other normative documents:

**EN 61000-6-3:2007+A1:2011**

**EN 61000-6-1:2007**

**Manufacturer Name & Address:**

Beyon Industry CO.,LTD

ADDRESS: No:3, Zhanhong Street, Dongcheng Road Dongshen town, Zhongshan City, China

TEL NO..+86-769-33238965


FAX NO.+86-769-89032220

**Date:**18.12.2017

**Name:** Richard Tang

**Title:** Sales Manager

**Signature:**
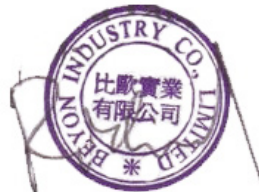
## Beyon Industry CO.,LTD

ADRES: No:3, Zhanhong Street, Dongcheng Road Dongshen town, Zhongshan City, China

TEL NO..+86-769-33238965

FAX NO.+86-769-89032220

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Biz, Beyon Industry CO.,LTD No:3, Zhanhong Street, Dongcheng Road Dongshen town, Zhongshan City, China

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